

MATTRESS

- ▶ **21 removable and independent Polyuréthane Ether cells** of 20 cm therapeutic AIR
- ▶ **4 therapeutic zones**
 - 2 static air cells located in head area
 - 4 cervico-dorsal cells
 - 10 air cells around the sacrum area with a static air base to ensure stability and avoid bottoming out
 - 5 independent heel relief air cells
- ▶ **Firm edge** for good stability when lifting and transferring
- ▶ **Side CPR** for deflation in less than 15 seconds
- ▶ **Transport plug:** System for opening and sharing 2 air channels (autonomy of 8 hours)



Heel relief



Side CPR


COVER SPECIFICATION

- ▶ Upper side: stretchable, breathable, impermeable to liquids and permeable to vapour
- ▶ Lower side: non-slip base
- ▶ **PROMUST PU HD cover**



PUMP

- ▶ Patented system of automatic and continuous calculation of inflation pressure according to patient's morphology
 - ▶ 3 modes: dynamic, static low pressure, care
 - ▶ Comfort settings: customization of comfort and seating position
 - ▶ Visual and audible alarms.
 - ▶ Indication for servicing
 - ▶ Serial number on the pump's back, simplified instruction and QR code on both sides.

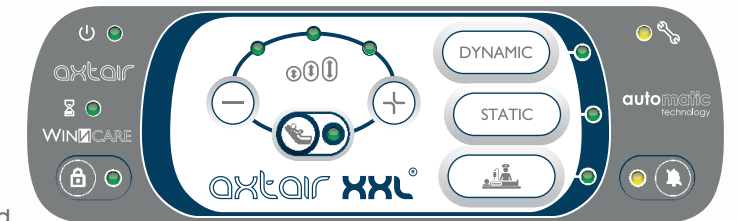

INDICATIONS

HELP IN PREVENTING PRESSURE ULCERS FOR PATIENT AT MEDIUM TO HIGH RISK OF PRESSURE ULCER DEVELOPMENT, BEDRIDDEN MORE THAN 15 HOURS A DAY.
 HELP IN TREATING EXISTING PRESSURE ULCERS FROM STAGE 1 TO 2 COMBINED WITH TECHNICAL AID DEVICES.
 PATIENT WEIGHT: 100 TO 270 KG/ 15.7 TO 42.5 ST

KEY BENEFITS

- ▶ **Effectiveness and therapeutic safety** for the patient
- ▶ **Easy to use** for caregivers
- ▶ Improved **quality of treatment and patient's quality of life**
- ▶ Firm edge for **good stability** when lifting and transferring
- ▶ **2 width** : 120 cm or 140 cm
- ▶ **Technical autodiagnostic** by QR code
- ▶ **Ready to use** in a carrying bag
- ▶ **Environmentally-friendly**
- ▶ **Eco-designed and manufactured in France**

DESCRIPTION	REFERENCE	DIMENSIONS / WEIGHT	PATIENT WEIGHT	FIRE STANDARDS	KEY FEATURES
AXTAIR XXL [®] (pump + overlay mattress) Promust PU HD cover For 120 cm wide bed	VAXT/XXL-120	195 x 117 x 20.5 cm	100 to 270 Kg	EN 597-1 & 2 GPEM D1 90 Classe D	Overlay mattress Side CPR 21 removable and independent Polyurethane Ether (300 u) cells of 20 cm therapeutic AIR 4 therapeutic areas : 2 static head celled, 4 cervico-dorsal cells, 10 cells around the sacrum area with air base, 5 heel relief cells Foam wedge Pump Dimensions: 22 x 25 x 11.5 cm Weight: 2.7 kg Technology based on alternating pressure 1 cell in 2 over a complete cycle from 9 to 14 minutes Patented system of automatic and continuous calculation of inflation pressure Comfort adjustment Visual and audible alarms Flow > 7 l / min Electrical power supply: 220-240 V - 50 Hz Class II type BF Average consumption: 5 Wh Pump decibels: Acoustic pressure according to NF EN3744: <35 dBA Electric and electromagnetic compliance according to the standards IEC 60601-1 Ed.3; IEC 60601-1-2
AXTAIR XXL [®] (pump + overlay mattress) Promust PU HD cover For 140 cm wide bed	VAXT/XXL-140	195 x 137 x 20.5 cm			2 years warranty

Cleaning instructions: Promust PU HD cover



* Behavior of the resistance of the compressor to disinfection with a solution of active chlorine at 1% (10,000 ppm) discontinuous use.

Photos not contractually binding

FP000102/03-2019

