

Luna

Ceiling hoist: Multiple solutions

The Luna ceiling hoist with its unique spherical shape and excellent functional features, is a true example of innovative Danish design.

Luna can be described as a fixed and portable ceiling hoist in one and the same unit, a truly unique feature. Changing Luna from fixed to portable mode, and vice versa, can be done within seconds by simply detaching and reattaching the unit to the track. No tools are needed. This unique two-in-one feature makes Luna the preferred ceiling hoist for care facilities as it provides a flexibility of choice, allowing the final decision of either fixed or portable unit to be made at a later point. This flexibility is a cost effective advantage for the institution.

Combine Luna with Ergolet's extensive programme of ceiling and wall mounted tracks, and you have a ceiling hoist system providing complete solutions to meet the latest requirements in the care sector.

Other features include Ergolet's patented spreader bar and sling hook system ensuring safe and simple sling attachment. And our comfortable ergonomic hand control, with its clear function buttons, is designed to help make patient transfer a safe and simple task for the carer.

Key standard functions

Class leading state-of-the-art design and functionality. Weight capacity of 200kg and 275kg.

Patented Ergolet sling hooks - simple and safe to use. Fixed or portable in a modular design.

Charging through the ergonomic hand control and wall/table charger.

Emergency lowering function.

Accessories

Complete selection of 2, 3 and 4 point spreader bars to suit standard and bariatric needs. Complete range of slings Weigh scales







Sling attachement: a minimum of downward pressure locks the sling hook. Pull upwards to unlock sling hook.



Mounts horizontally in a fixed or portable mode without using tools.

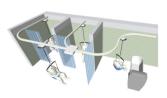



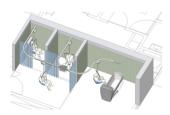

Room to room transfer.



Version	ltem-number
Luna 200 kg (white/grey)	015-01910 / 015-01912
Luna 275kg (white/grey)	015-01936 / 015-01937
Luna room to room 200kg (white)	015-01922
Luna room to room 200kg (grey)	015-01920

TECHNICAL SPECIFICATION

Max. capacity 200 kg / 275 kg Length of lifting strap 2200 mm

Weight 6,7 kg w/o spreader bar

Measurements 320 x 130 mm

Battery operation 36 V, lead acid batteries Charging 100 - 240 V, 0,5A.

Chassis Aluminium

Housing High impact, moulded ABS

Noise Lwa during lifting/lowering without load is 49dB(A)

Lwa during lifting/lowering at maximum load is 53dB(A)

Lifting speed down

Lifting speed up

30 mm/sec with 89 kg load

No. of lifts per charge


44 mm/sec with 89 kg load


50 lifts of 0,5 m with 100 kg

ISO ISO 10535

- Soft start and stop motor control
- Electrical emergency lowering
- Electrical emergency stop
- Low battery indicator (visual LED)
- Charging indicators
- Power on indicators
- Overload protection
- Emergency brake system

Fixed Portable

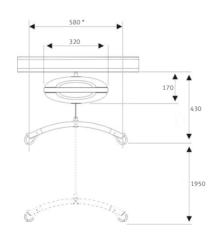

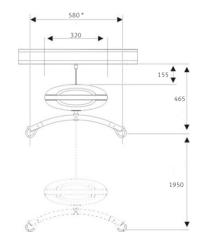

 $^{^{\}star}$ The size of the spreader bar varies depending on the model.

Consult the Technical manual for further specifications. All figures and measurements are for guidance only and may differ in actual situations. Ergolet reserves the right to changes without notice.

Product photos may include accessories. Consult your dealer.

Ergolet
Taarnborgvej 12C
4220 Korsoer
Denmark
tel. +45 70 27 37 20
e-mail: ergolet@ergolet.dk

