



Therapy & comfort

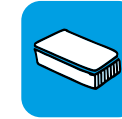

Score on the Braden scale

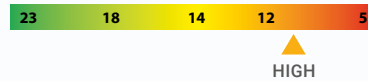

axtar automorpho<sup>®</sup>  
plus

## THERAPEUTIC INDICATIONS

Helps prevent pressure ulcers in **MEDIUM** to **HIGH RISK** PATIENTS.  
Patients confined to bed for over 15 hours a day.  
Helps treat **stage 1 to 4** pressure ulcers used in conjunction with a local pressure relief or positioning aid.

30 ≤ 180 kg



FIRE REGULATIONS  
EN 597 parts 1 & 2  
GPEM/CP D1-90 Class D



3 bed widths available 90 | 100 | 120 cm


## KEY CHARACTERISTICS

- ▶ Patented system for automatic and continuous adjustment of inflation pressure depending on patient morphology
- ▶ Built-in heel rest
- ▶ Comfort settings for all modes and patient positions
- ▶ 3 widths available: 90, 100 and 120 cm
- ▶ QR code for quick and easy access to product information
- ▶ Delivered in transit bag - remove and it's ready to use
- ▶ Environmentally friendly design
- ▶ 2 covers available: Promust PU HD and Promust CIC
- ▶ Mattress weight: ~ 10 kg (+ 2.7 kg compressor)
- ▶ Integrated pass cable management
- ▶ Rewash and barrier machine compatible



## KEY BENEFITS

- ▶ Improved patient care and quality of life
- ▶ Easy to use for care personnel

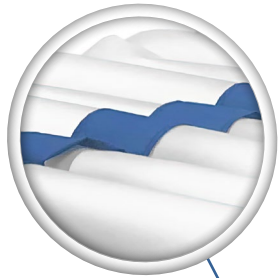


## COMPARTMENTS

- ▶ 18 polyurethane compartments height 12 cm thickness 300  $\mu$ , filled with therapeutic air



## CPR (Cardio Pulmonary Rescue) VALVE

- ▶ Quick and easy opening
- ▶ Mattress deflates in under 15 seconds in emergencies\*  
*\* With a patient weighing 80 kg and lying flat*

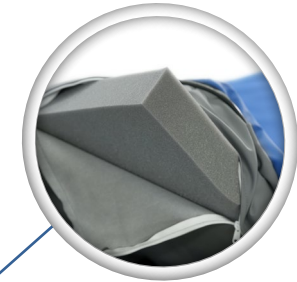



## FOAM BASE

- ▶ Helps keep the patient safe in the event of undesirable incidents




## HEEL REST

- ▶ Déflation of one or more compartments at foot of mattress relieves pressure on the heels

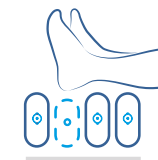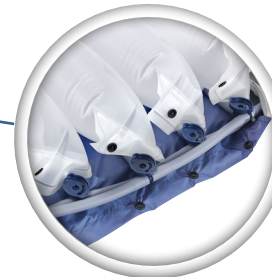


Compartment is deflated then withdrawn (remove the cap)

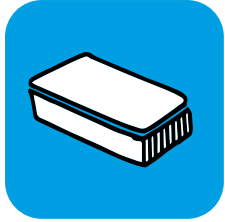

# axtair automorpho<sup>®</sup> plus

## COMPRESSOR

- Patented system for automatic and continuous adjustment of inflation pressure depending on patient morphology
- 3 operating modes: dynamic, slow pressure, care
- Management of seated position
- Comfort adjustment in all modes and patient positions
- Visual and audible alarms
- Serial number on back of compressor, simplified notice
- QR code 'SCAN ME!' for quick and easy access to information
- Maintenance management with the AIRCARE program

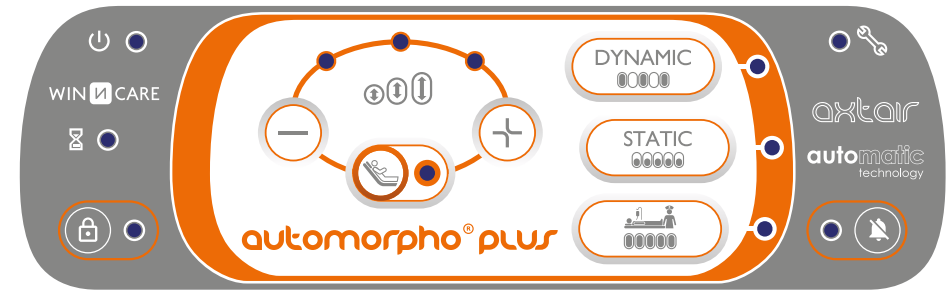

## TECHNICAL CHARACTERISTICS OF THE COMPRESSOR


DIMENSIONS	22 x 25 x 11.5 cm
DYNAMIC MODE (ALTERNATING)	1 compartment in 2
TIME TO COMPLETE FULL CYCLE	9 to 14 minutes
PATENTED TECHNOLOGY	Patented automatic control system
ADJUSTABLE SETTINGS	Personalised comfort
INFLATION SPEED	> 7 l / min
POWER SUPPLY	220-240V - 50Hz
CLASS	Class II type BF
AVERAGE POWER CONSUMPTION	5Wh
SOUND PRESSURE UNDER NF EN3744	< 35 dBA IP41
ELECTRIC AND ELECTROMAGNETIC COMPATIBILITY	IEC 60601-1 Ed. 3; IEC 60601-1-2
PROTECTION INDEX	IP42





## MATTRESS-COMPRESSOR ELECTRICAL INTERFACE

- ▶ Equipped with a transit plug for combination of two channels in transit mode for patient comfort and 8 hours' autonomy

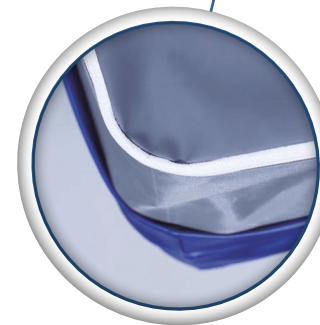



## BOTTOM COVER

- ▶ Fitted with a non-slip base

## COVER

- ▶ Upper surface: bi-elastic, impermeable to fluids and permeable to water vapour
- ▶ Lower surface: anti-slip base
- ▶ 2 Covers available:  
**PROMUST PU HD**  
**PROMUST CIC**



COVER	WIDTH OF BED (CM)	INDIVIDUAL CARDBOARD BOX
	90	VAXT4/AUTO-P
	100	VAXT4/AUTO100-P
	120	VAXT4/AUTO120-P
	90	VAXT4/CIC-P
	100	VAXT4/CIC100-P
	120	VAXT4/CIC120-P